

IM-5398H

Instruction

CONTROL PANEL





Dear customer,

Thanks for your choosing our products.

For the best capability and safety, please read the manual before operation.

Please reserve the manual for your reference.

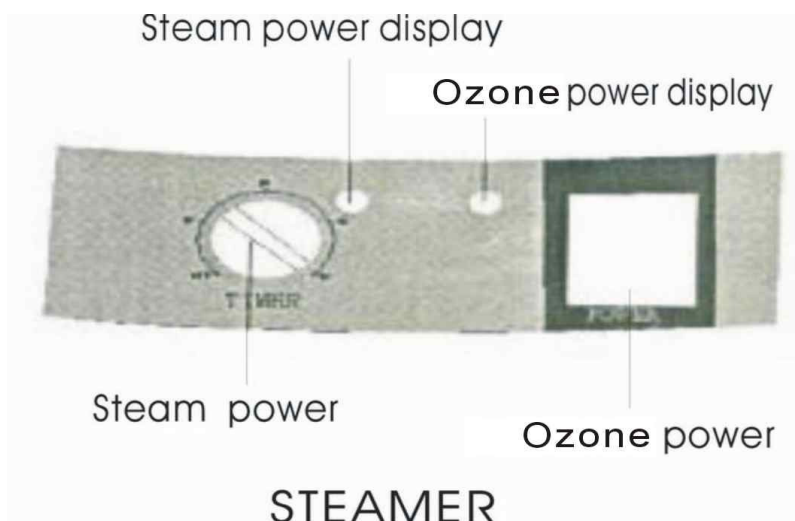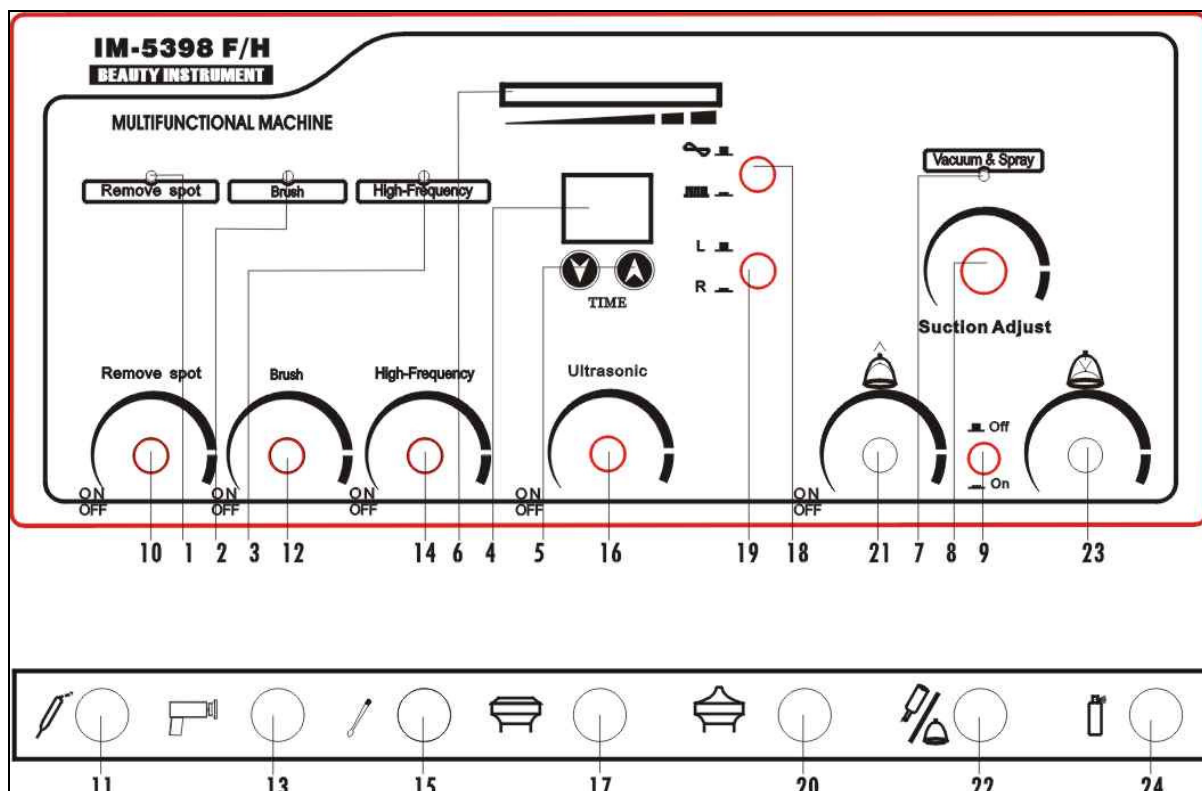

- | | |
|--|-------------------------------------|
| 1. Spot Remover Display | 2. Brush Display |
| 3. High-Frequency Display | 4. Ultrasonic Time Display |
| 5. Ultrasonic Time Adjuster | 6. Ultrasonic Intensity Display |
| 7. Sprayer & Vacuum Display | 8. Suction Adjust |
| 9. Switch for Sprayer& Vacuum | 10.Remove Spot Intensity Adjust |
| 11. Remove Spot Output | 12. Brush Intensity Adjust |
| 13.Brush Output | 14. High Frequency Intensity Adjust |
| 15. High Frequency Output | 16. Ultrasonic Intensity Adjust |
| 17. Big Probes Output | 18. Ultrasonic Output Wave Switch |
| 19. Big / small Probes Output Switch | 20. Small Probes Output |
| 21. Patter Breast Care Switch/ Suction Time Adjust | |
| 22. Vacuum/Patter Breast Care Output | |
| 23..Patter Breast Care Relaxing Time Adjust | 24. Spray Output |

HIGH FREQUENCY SPOT REMOVER

The High frequency Spot remover is through the high frequency electricity rate coagulation, urges the epidermal cell necrosis, which can in instantaneous thoroughly eliminate all kinds of spot both on the face and on skin surface, such as old age spots, freckles, granules, tattoos, moles, etc. The effect is obvious with no bleeding or scar left.

Procedure:

1. Sterilize the devices with 75% alcohol.
2. Plug in the power supply, and plug the spot remover into relevant slot.
3. Switch on the power, and adjust the density to suitable degree.
4. Grasp the spot remover and carefully point on the spotted mole, until seeing Melanin excreted as waste.
5. Apply medicine for disinfections on the treated surface.
6. Please grease antiphlogistin after the treatment.

Cautions:

1. Before using, the spot remover and the surface should be cleaned and sterilized.
2. After every usage, remove the power supply, and the spot remover should be clean and sterilized
3. When operating, do not stab the spot remover in the skin too much, but only for recommended duration to avoid burns on the dermis.
4. If spot can not be removed during one usage, do not use again until the scab falls off.
5. Removing no more than 5 spots is recommended in one usage.
6. After the usage, it is recommended the client not to wet the operated surface within 22 hours, do not eat ginger and anything with dark color.
7. Do not tear away the scab, must let it fall off to avoid infections and scars.
8. Do not expose intensely under the sunshine.

OZONE HOT STEAMER

Function:

- ⊙ Accelerate blood circulation and support cell water.
- ⊙ Expand pore and deeply clean.
- ⊙ Remove aging horniness and grime.
- ⊙ Diminish inflammation and sterilize wound.
- ⊙ Improve skin metabolism and rejuvenation epidermis cell.

Procedure:

1. Use measuring cup to fill distilled water in the affusion gate and make sure the water should be higher than the affusion line.
2. Connect the wire and check the power.
3. Press the red key, after 5 to 6 minutes, the steam will come out. This is the normal steam; and then press the green key, the steam together with the ozone will come out, which have sterilization function.
4. Please cover eyes with wet cloth patch and don't face the steam directly before it becomes steadily.
5. Adjust the distances between the spray spout and client's face according to the nature of the skin. The steam should come through the face from upside to underside.
6. Please turn off the machine and power when finished.

Time and Distance:

The treatment time and distance between spray spout and client's face should be adjusted according to the nature of skin.

- (1). Normal Skin: distance----25cm~30cm; time (ozone steam)----3~5 minutes
- (2). Greasy Skin: distance----20cm~25cm; time (ozone steam)----5~8 minutes
- (3). Sore Skin: distance----20cm~25cm; time (ozone steam)----8~10 minutes
- (4). Dry Skin: distance----30cm~35cm; time (ozone steam)----3 minutes
- (5). Sensitive Skin: distance----35cm; time (normal steam)----5~8 minutes
- (6). Spot Skin: distance----30cm~35cm; time (normal steam)----5~8 minutes

minutes

(7). Telangiectasis Skin: distance----20cm~25cm; time (normal steam)----5~8 minutes

For all the treatment with ozone, the total time should be within 10 to 15 minutes. Please don't do it over 15 minutes.

Attention:

1. Avoid spray directly to the nostril, as it is easy to cause breath problem or uncomfortable feeling.
2. The treatment time should be less than 15 minutes.
3. If spray arm is loose, can firm it as the picture instruction
4. Water cup should be cleaned periodic.
5. When adding water, please turn the power off first.
6. When finish, please turn off the power.
7. If furring occurs to clog the steam, please pour in some vinegar and keep for about 15 minutes.



SKIN VACCUMER & SPRAYER

1. Basic principle of Vacuum & Sprayer

When air pump is working, we can use vacuum to suck out the dirty things in the skin. The air pump make strong air circulation, connect the spray bottle to the vent hole, then can spray cosmetic water. We can adjust the vacuum intensity according to specific situation.

2. Function of Vacuum & Sprayer

- (1).Remove dirty things in the skin, give the skin a fresh breath. Stimulate the skin organization, strengthen the tenderness of skin.
- (2).Accelerate the blood circulation, lead blood to scarfskin, supply nutrition to the skin, promote lymphatic circulation, drain out toxin inside the skin.
- (3).The spray function can evenly spray cosmetic water on the skin and

make the skin absorb it very well, to help shrink the pore.

3. Main tools

Alcohol cotton, cotton stick, sterile cotton, cosmetic water, vacuum instrument and its accessory

4. Operation procedure

Vacuum

- (1). Connect the sterilized glass tube to the pipe which connected with vacuum output.
- (2) Start the machine and feel the suction intensity on the back of hand, adjust the intensity according to specific demand.
- (3). There are three ways to make vacuum: discontinuous vacuum, continuous vacuum, strong vacuum. When make treatment to the oily skin, hard and skin, need use much stronger vacuum. When make treatment to xeransis skin, sensitive skin, senile skin, need little lighter vacuum.
 - a. Discontinuous vacuum. During operation, the middle finger should cover and release on the small vent hole of glass tube in a frequent and rhythm way. The vacuum should be not too strong. This kind of operation is better to treat tender skin, loose skin and thin skin.
 - b. Continuous vacuum. The middle finger cover the small vent hole of glass tube, and move together with the tube for a while then release the finger. The vacuum need much stronger. This kind of operation is better to treat oily skin, hard and thick skin.
 - c. Strong vacuum. The middle finger cover the vent hole of glass tube and don't release, put the tube mouth on the skin, suck and release on the skin frequently. The vacuum intensity need strong. This kind of operation is better to treat oily skin, such parts as nose tip, nose side and the parts with comedones.

Spray

- (1) Connect the spray pot that filled with nutrient fluid to the pipe which connected with spray output.
- (2) Start the machine, put a figure to cover on the round hole of the spray pot top, check if there is fluid spray out.
- (3) Until your client close his eyes, please spray the fluid with the direction

from forehead to mandible and keep a suitable distance with the face.

5. Operation Notice

- (1) Should not move the vacuum tube too slow, should not vacuum at a part for too long time
- (2) According to client's skin nature to adjust the suction intensity, should try it on the hand back at first, then apply to the face part.
- (3) Should match with the steamer to make the treatment.
- (4) Around the eyes, the skin is too thin, should not apply the vacuum.
- (5) Before spray, please ask the client close the eyes, and the spray direction is from forehead to mandible and keep a suitable distance with the face.
- (6) A glass without the small vent hole, we can use it to stimulate lymph, but should with lighter vacuum.
- (7) All the glass should be sterilized after finishing operation.
- (8) When treating the nose comedones, we should put some cotton into the glass in order to avoid the waste is suck into the machine which is harm for it.
- (9) Please pay attention to the output of vacuum and spray, make sure with the right output when you use vacuum or spray. If not, it will lead to breakdown.

HIGH FREQUENCY ELECTROTHERAPY

The high frequency electrotherapy is an important beauty instrument. It is through the glass electrode tube to produce high frequency current, and apply on the face, thus achieve the treatment result.

1. It helps to calm the nerve and has analgesia function.
2. It produces the ozone for sterilization, can help wound to heal and to restrain pores.
3. It can promote skin's secretion and improve the PH value, which makes the skin neutrality.
4. It can accelerate blood circulation and improve metabolism.

5. It has the function of anti-aging, anti-wrinkle and tightening.

Operation Procedures:

1. Indirect Method:

- (1). Please screw in the glass tube into high-frequency handle..
 - (2). Please ask the client to hold the glass tube after the hands smearing some talcum powder. Turn on the machine and the electricity will get through from hands to the whole body.
 - (3). Beauticians do massage on client's face from neck to chin and from visor to forehead.
 - (4). Please adjust the intensity to "0". Make sure that the machine is turned off and puck the glass tube out of the handle.
- This method is fit for dry nature skin and aged skin.

2. Direct Method:

- (1). Please grease essence and cream on the face first and then the beautician hold the metal bar to put the glass tube on client's face. Turn on the machine.
 - (2). Massage with the glass tube like "Z" or helix on client's face, do the treatment as follows, forehead→nose→right side of the face→chin→left side of the face→nose→forehead
 - (3). Please adjust the intensity to "0". Make sure that the machine is turned off and puck the glass tube out of the hand piece.
- This method is fit for greasy nature skin.

3. Sparkle Method:

- (1). Cover client's eyes with a piece of the wet cloth patch and touch the face with the glass tube on the part of inflammation. For one time, the treatment should be less than 10 seconds.
 - (2). It's normal if there are sparkles come out while the glass tube touches the face.
 - (3). Please adjust the intensity to "0". Make sure that the machine is turned off and puck the glass tube out of the hand piece.
- This method is fit skin with sore, inflammation and wounds.

Cautions:

1. To tender skin, please cover a piece of voile to avoid the excess stimulation.
2. Do not place on one part over 10 seconds.
3. Pregnant woman and people with pacemaker are unfit for this treatment.
4. Please make sure that install the glass tube into the high-frequency handle before turning on the machine. Adjust the intensity to “0” after finish.
5. Skin with freckle and splash is unfit for this treatment.

BRUSH**Working theory**

Making use of the rotary brush with high speed for scrub both on face and body to improve the skin by removing spot, wrinkles and aged surface. Because it's different from the medical treatment, as this rotary brush with accessories can offer you for gentle exfoliation and delicate stimulation in a safer way, which will give your client a radiant complexion with no harm to the skin.

Product Character

1. No side-effect and won't leave scar and spot on the surface
2. Do not affect the daily work and social life
3. Can be used on all kinds of skin

Operation Instruction:

4. Connect the brush and the handgrip to the mainframe
5. According to the skin nature, please choose different brushes
6. Turn on the working power, adjust the proper frequency
7. The proper working time for the whole is settled in 30 minutes
8. Turn off the power after finished and sterilize the brush

Attention:

In one to two days after the treatment, you should pay attention to the tips as follows:

1. Do not do facial or shave
2. Do not go to the hot and public place after the treatment, such as Sauna, hot smoke, etc.

ULTRASONIC

The ultrasonic function is suitable for any skin problem nature, with a simple and easy operation; it is the common instrument among esthetical professionals.

The ultrasonic acts on the human body and is absorbed by body, then transfers the sound energy into heat energy which provides stimulation for human skin, accelerates the metabolism, circulates the blood and lymphatic, activates the deep inside cell and improves the defense system.

Functions:

1. The ultrasonic helps to discoloration, and reduces the damaged area.
2. With its massage effects, it helps the skin to better absorb nutrition.
3. It accelerates the partial blood circulation, improves the metabolism.
4. It dissolves hypodermic fat, accelerates the absorption, cause to discompose the agglomerated excessive water and fat, which can improves eye bag and black eye circle problems.
5. It heals acne skin.
6. The ultrasonic wave produces the temperature effect, and helps to calm the nerve and has analgesia function.

Operation Procedures:

1. Big probe for face, continuous wave, for 15 minutes.
2. Small probe for eyes, pulse wave, 5 minutes for each.
3. Duration, intensity and wave type should all depend on skin nature.
 - a. Extra dry skin: mature, spot and aged skin

- Duration: 13-15 minutes; Wave: Continuous
- b. Dry skin: Thinner, red capillary skin
Duration: 10-12 minutes; Wave: Continuous
 - c. Normal skin: Healthy skin
Duration: 8 minutes; Wave: Continuous
 - d. Oily skin
Duration: 8 minutes; Wave: Continuous
 - e. Acne skin
Duration: 8-10 minutes; Wave: Pulse
 - f. Sensitive skin
Duration: 6-8 minutes; Wave: Pulse
 - g. Eye (each)
Duration: 4-5 minutes; Wave: Pulse

Cautions:

1. Do not touch the eyeball with the probe.
2. Apply sufficient essence on the surface before starting the treatment.
3. Do not point and stay at any part over 15 minutes.
4. Clean and sterilize the probe, should put in cool, dry place to prevent oxidation
5. Disconnect the power when finish.
6. Use pulse wave for inflamed and damaged skin.

MAGNIFYING LAMP

With easy to move and fix position features, the Magnifying Lamp includes a superior refractive magnifier and a built-in illumination system that allows you to examine skin for great detail.

Operations:

1. Plug in, and put through electrical outlet;
2. Open the on-off knob and let the magnifying glass point to the skin which need for checking and operating on;
3. Adjust the distance between Magnifying Lamp head and skin.

Then you shall fix the knob. (Annotate: you'd better use some body shaping product on part of body where need examined);

4. After operation, turn off the on-off knob and close up electrical outlet;

Attentions:

1. If the power electrical wire was damaged, to avoid the danger, please contact with the manufacturer or after-sales services suppliers or other professionals to help you replace it.
2. Do not lay it in the place of dank or moisten environment, or else it will bring short circuit and fire etc.
3. Don't take it down when the apparatus is electrify.
4. You can not touch it if it was dipped by liquid;

WOODS LAMP

An advanced way to analyze skin is by using a skin scope. Backed by the technology of facial topography, this skin scope can analyze skin condition and detect and interpret skin problems. It has the ability to identify skin types that co-exists on an individual's face. Skin diseases and problems can be easily detected using a skin scope. Uses black light to illuminate imperfections and dryness in the client's skin. This machine is able to examine skin in minute detail. Indications of Color:

Color and Skin:

2. White with cyan: Healthy skin
3. Yellow with cyan: Greasy skin
4. Purple with cyan: Dry skin
5. Dark purple: Highly dry skin, Dehydrated Skin
6. Orange: Acne greasy skin
7. Light yellow: Acne fester skin
8. Brown and dark brown: Pigmentation skin
9. Purple: Sensitive skin
10. White fluorescence: Thick Corneum Skin

Operation:

1. Please cover eyes with wet cloth patch after skin clean
2. The beautician holds the lamp and take the lamp towards client's face. The distance between face and lamp should be 15cm~20cm.
3. Looking the color of the skin reflect on the lamp
4. After finish the checking, please put away the lamp

Caution:

1. Please cover client's eyes with web cloth patch before start.
2. The testing time should be less than 2 minutes.
3. Please make sure that the distance between face and lamp should be no less than 15cm.

PATTER BREAST CARE

Breast is valued for all women and keeping it healthy and plump is the dream of all ladies. This product is using excise way to strengthen breast muscle and active tissue, thus to achieve the effect of breast plumping and care.

Operations:

1. Turn on suction switch and the machine starts sucking but no releasing.
Using this function to remove black heads、do medical jar extracting treatment etc. Connecting the glass vessel to the sucking conduit and press the hole on the glass vessel to produce interval suction result.
2. Tie the cups to the three-tie glue pipe and set at maximum, fix cups on the breast for the clients (with some gel ointment or flavor oil)and inquire client's feeling to adjust frequency, intensity and time(normally set at 30 minutes, if necessary, could prolong it) to achieve eligible result.
3. Clean the cups and glass vessels and put away after usage.

Important Safety Notice

1. Please don't disassemble the equipment or try to do some other operation, which have not been indicated in our instruction. All the repairing work should be done by our admitted professional personnel
2. Please don't assemble or operate this equipment while you are near the water or your hands are wet. Notice that don't pour the liquid on it.
3. Please pull out the plug from the socket when you met such conditions as below. And keep in touch with professional maintenance Station.
 1. The equipment touched liquid
 2. The equipment has abnormal smell, smoking or loud and strong noise.
 3. Cable is broken
 4. The equipment is fallen down and broken
4. Please don't put anything on the cables, don't place the equipment at any place where someone can step on the cable.
4. Please pull out the power cable after using it for your safety
6. Don't put anything into the equipment's jack, otherwise it may cause fire or touch the electricity. If something enters the equipment, please don't pick it out by yourself, contract distributor or our company immediately
7. Please don't put the equipment on the anything that can't stand steadily, avoid it turns over and damage.
8. Our company reserves the right to modify the information in this introduction before notify! And keep the explanations right for the above.

Technical Standard :

Characteristics

·Power Supply	·□220V/50Hz . □110V/60Hz □240V/50Hz
· Power	·880W
·Operation Temperature	·5℃-40℃
· Humidity	·80% without condensation
·Dimension	·103cm (L) ×68cm (W) ×41cm (H)
·Weight	34Kg (including accessories)

List of accessories

Description	Unit	Qty	Remark
Mainframe	Set	1	
Steamer (only head)	Set	1	
Magnifying Lamp	Set	1	
Wood Lamp	Set	1	
Vacuum Glass	Set	1	With 3 pieces
Brush Probe	Set	1	
Brushes	Set	1	With 5 pieces
Remove spot handle	Pc	1	
Spanner	Pc	1	
Ultrasonic big probe	Pc	1	
Ultrasonic small probe	Pc	1	
High-frequency handle	Pc	1	
High-frequency glass tube	set	1	With 3 pieces
Spray pot	pc	2	
Breast care plastic cup	Set	2	70#/90#
Trifurcate Plastic Pipe	Strip	1	Breast care
Salver	Pc	1	
Spray & Vacuum connect pipe	Strip	1	
Filter	Pc	1	
Fuse	Pc	2	
Power Wire	Strip	1	
Instruction	Pc	1	